



AIR-TRAQ



Escambia County, Florida - ALS/BLS Medical Protocol

The AIRTRAQ is a disposable optical laryngoscope visual guiding system for routine and complex airway use.

INDICATIONS

Adjunct for Endotracheal Intubation

CONTRAINDICATIONS

The AIRTRAQ should not be used w/< 7mm or >8.5mm ET Tubes.

INSERTION STEPS

1. Turn the light on 30 seconds prior to use to allow lens to warm for de-fogging.
2. Lubricate and install the endotracheal tube in the lateral channel.
3. Insert the AIRTRAQ into the mouth, sliding it along the midline.
4. Look through the eyepiece as insertion is continued, to identify airway structures.
5. Identify the epiglottis and place the tip of the AIRTRAQ into the vallecula.
6. Gently lift and maneuver the AIRTRAQ to center the vocal cords into view.
7. Advance the endotracheal tube, watching the tip and cuff pass through the vocal cords.
8. Inflate cuff prior to removal of Endotracheal Tube from Airtraq .



Approved: 10/01/2011
Version: 3.1

Page 1 of 4

Approved by:


Charles Neal, D.O. Medical Director



AIR-TRAQ



Escambia County, Florida - ALS/BLS Medical Protocol

REMOVAL STEPS

After endotracheal tube cuff inflation, separate the AIRTRAQ from the endotracheal tube and remove it while holding the endotracheal tube in place.



Approved: 10/01/2011
Version: 3.1

Page 2 of 4

Approved by:


Charles Neal, D.O. Medical Director




AIR-TRAQ

Escambia County, Florida - ALS/BLS Medical Protocol

AIRTRAQ®
OPTICAL LARYNGOSCOPE

**Refer to full IFU for detailed information and user tips.*

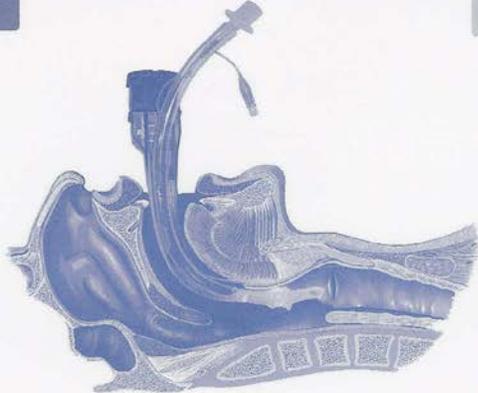


Anatomically-shaped laryngoscope with two separate channels. The **Optical Channel** contains a high definition optical system. The **Guiding Channel** holds the endotracheal tube and directs it towards the vocal cords.

Built-in **Anti-Fog System** and light source.

AIRTRAQ			
SIZE	PART NO	COLOR	ETT SIZE
Regular	ATQ-011	Blue	For use with ETT 7.0 – 8.5 Minimum patient mouth opening: 18 mm
Small	ATQ-021	Green	For use with ETT 6.0 – 7.5 Minimum patient mouth opening: 16 mm

Clip-On Video System Accessory
Part No ATQ-032

www.kingsystems.com
800-642-5464 or 317-776-6823


AIRTRAQ is a registered trademark of Prodel Meditec S.A. U.S. Patent No 6,843,769. PIFM-05, IIT-05 05/06



Approved: 10/01/2011
Version: 3.1

Page 3 of 4

Approved by:

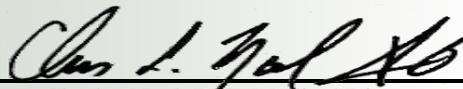
Charles Neal, D.O. Medical Director



AIR-TRAQ



Escambia County, Florida - ALS/BLS Medical Protocol

AIRTRAQ[®]
OPTICAL LARYNGOSCOPE
*Refer to full IFU for detailed information and user tips.



insertion steps

- Turn the light on.
- Lubricate and install the endotracheal tube in the lateral channel.
- Insert the AIRTRAQ into the mouth, sliding it along the midline (**Fig 1**).
- Look through the eyepiece as insertion is continued, to identify airway structures (**Fig 2**).
- Identify the epiglottis and place the tip of the AIRTRAQ into the vallecula.
- Gently lift and maneuver the AIRTRAQ to center the vocal cords into view (**Fig 3**).
- Advance the endotracheal tube, watching the tip and cuff pass through the vocal cords (**Fig 4**).

removal steps

- After endotracheal tube cuff inflation, separate the AIRTRAQ from the endotracheal tube and remove it while holding the endotracheal tube in place (**Fig 5**).

- OVER -

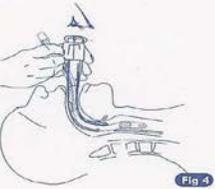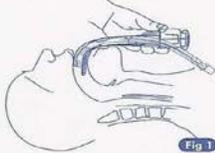




Approved: 10/01/2011
Version: 3.1

Page 4 of 4

Approved by:

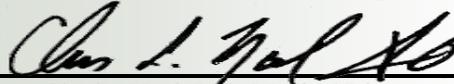
Charles Neal, D.O. Medical Director